


EXHIBIT B

APL RESEARCH CENTER
ANALYTICAL RESEARCH DEPARTMENT

SAMPLE ID : PR(094)98
FILE NAME : D:\DATA\████████
METHOD NAME : D:\METHODS\CEFDINIR_gradient.met
ACQUIRED TIME: ██████████

VIAL NO: 1
INJ. VOL : 20 μ l
INSTRUMENT ID: ARE_058
USER REF: kk -MD()

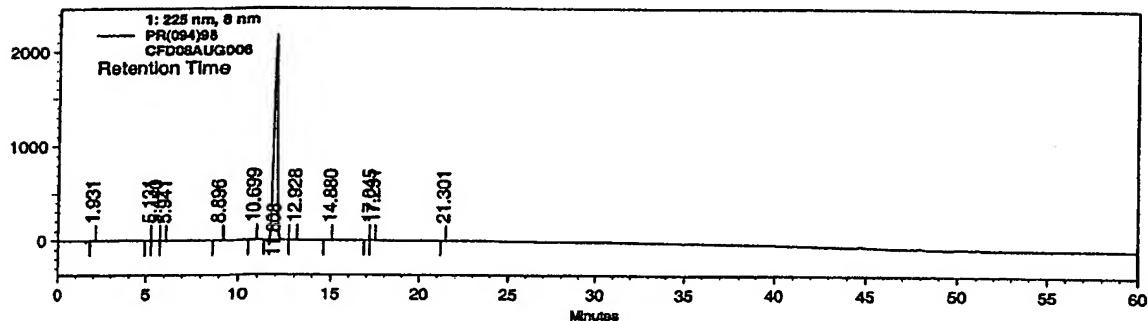

1: 225 nm, 8 nm	Retention Time	Area	Area Percent	Relative RT
	1.931	13912	0.05	0.16
	5.131	7636	0.03	0.43
	5.440	9639	0.03	0.46
	5.941	11407	0.04	0.50
	8.896	14601	0.05	0.75
	10.699	21918	0.07	0.91
	11.808	29201853	99.40	1.00
	12.928	10317	0.04	1.10
	14.880	11330	0.04	1.26
	17.291	21316	0.07	1.45
	21.301	39769	0.14	1.47
		13306	0.05	1.81
	Totals	29377004	100.00	

HPCC LABORATORY NOTE BOOK

12

OK

Product: Caffeine

Project: Caffeine

S/N: ARF-058

Method for: Chromatographic purity

Here Column:-MP: (0.02M NaHPO₄ + 0.4mM EDTA) pH 6 : (Acet + MeOH) 8:2

60	98	02
50	90	10
30	65	35
45	30	70
60	30	70
61	98	02
70	98	02

Column: HyperSil ~~box~~ C18 Crossed 5μ Flow: 1.0 ml/min

Det: 225nm Resolution: 60 mins T = 40°C

Observation & Conclusion:-

RT PR(94) 98 Sample Analyzed purity 99.4%

Caffeine: 11.9

Amide: 42.0

AVN: 4.5

All closely eluting impurities are well separated
 Method is good but baseline drift towards neg.
 is observed, it should be improved.

HPLC INSTRUMENT LOG BOOK

ၬ

CONTINUATION 17

1ST AVAILABLE COPY